

Brief Operating Instructions – Blood Glucose Measurement with BG★STAR®



- ★ **No coding** necessary
- ★ **Sample size** – 0.5 µL
- ★ **Test duration** – 6 seconds*
- ★ **Memory** – 1,865 measured values

* average length of test

These instructions are intended solely as a guide. Please read through the detailed operating instructions carefully before using the BGStar® for the first time.

1



Insert test strip

- ★ Insert the test strip into the test strip slot with the contacts showing upwards. The meter will switch on automatically.
- ★ **WARNING:** To ensure accurate results, please wash your hands with warm soapy water and dry them before each test.

2



Test site

- ★ Follow the instructions to prepare the lancing device. You will find the instructions for the lancing device in the BGStar® pack or overleaf.
- ★ Press the lancing device against the place where you would like to pierce the skin. Press the trigger button. To take blood from other parts of the body, please read the operating instructions.

3



Apply the test strip to the blood sample

- ★ The test strip draws the blood through the test area into the strip.

4



Remove the test strip from the blood sample when the bleeper sounds

- ★ Your meter will beep once if no further blood needs to be drawn into the sample area of the test strip. The progress bar and the "1-2-3" display indicate that the meter is calculating your results.


5




Final steps

- ★ Your test is completed once your blood glucose result, with the time and date, is shown on the display. By pressing the "Confirm" key, you confirm the mealtime marking and the marking will stop flashing.
- ★ Remove the used test strip. The blood glucose meter switches off automatically. Pull the cap off the lancing device and place the lancet cap on the lancet.
- ★ Remove the used lancet and dispose of it properly.


Preparing the lancing device

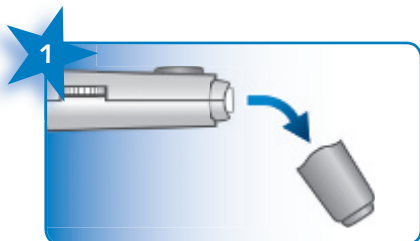


Pull the cap off the lancing device.

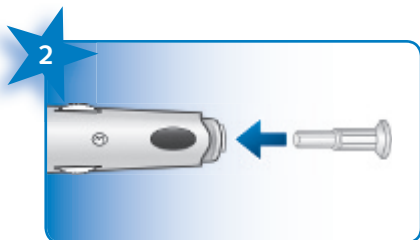

Place the new lancet firmly on the lancet holder.

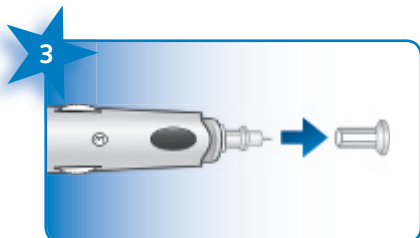


Twist the lancet cap off.

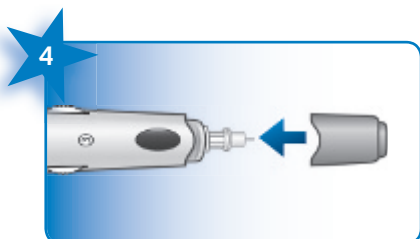

Place the cap back on the lancing device.



Adjust the piercing depth.

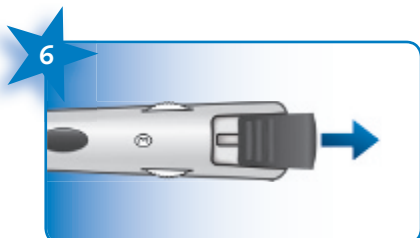

Pull on the tensioning device until it clicks.



www.bgstar.de

Your Diabetes Service Centre

Sanofi-Aventis Deutschland GmbH
ServiceCenter Diabetes
Industriepark Höchst
65926 Frankfurt am Main